

TMSi PC-Driver Installation Manual



Table of Contents

General information	4
How to use this manual	4
Minimum system requirements	5
Supported hardware	5
Supported operating systems	5
Additional information	5
Uninstall previous versions	7
Installation	7

General information

- This manual contains the information needed to install the TMSi PC-Driver v7.x.xxx.x.
- The TMSi PC-Driver is designed to be used with a PC with a Windows operating system.
- The screen shots in this manual are taken during the installation of the TMSi PC-Driver on a PC with Windows 8.

How to use this manual

- Read this manual carefully before installing and using the equipment.
- Complete the actions in the given sequence.
- Keep a copy of this manual near the equipment.

Minimum system requirements

Processor: Celeron 2.0GHz

Memory: 1GB RAM

Supported hardware

The silver-colored label on the product shows the serial number. This label can be located on any side of the product. The first four digits of the 10-digit serial number are the device code. Devices with the following device codes are supported:

- 0107, 0111, 0120, 0121, 0700 = Refa
- 0112, 0118, 0122, 0128 = QuickAmp
- 0207 = Porti
- 0208, 0934, 0938, 0928 = NeXus
- 0392 = Fusbi
- 0393 = Synfi
- 0710 = Mobita
- 0925, 0926, 0928, 0931 = Mobi(Mini)

Supported operating systems

The following operating systems are supported:

- Windows 7
- Windows 8.1 (except for Mobita)

Additional information

For information on changes between the TMSi PC-Driver and its predecessor WDM driver v6 refer to TMSiSDK.chm in C:\Program Files \TMSi PC-Driver\ExampleCode.

Uninstall previous versions

Warning



Before starting the un-installation of the TMSi PC-Driver make sure no TMSi devices are connected to the PC!

Windows 7 & 8: Uninstall previous version of the TMSi PC-Driver through the Control Panel → Programs → Programs and Features.

Installation

Warning



Before starting the installation of the TMSi PC-Driver make sure no TMSi devices are connected to the PC!


Warning




Uninstall previous versions of the TMSi PC-Driver before starting the installation of this TMSi PC-Driver !

The installation starts after double-clicking 'setup.exe':


Step 1: If User Account Control dialog pops up, click 'Yes'.

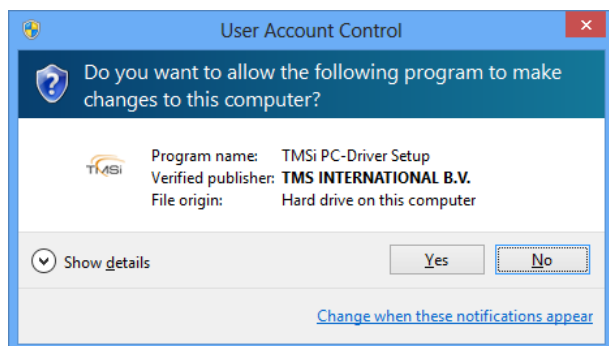


Step 2: Select Setup Language and click 'OK'.

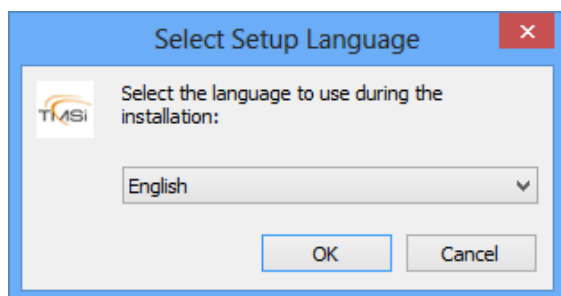


Step 3: Setup – TMSi PC-Driver, click ‘Next’.

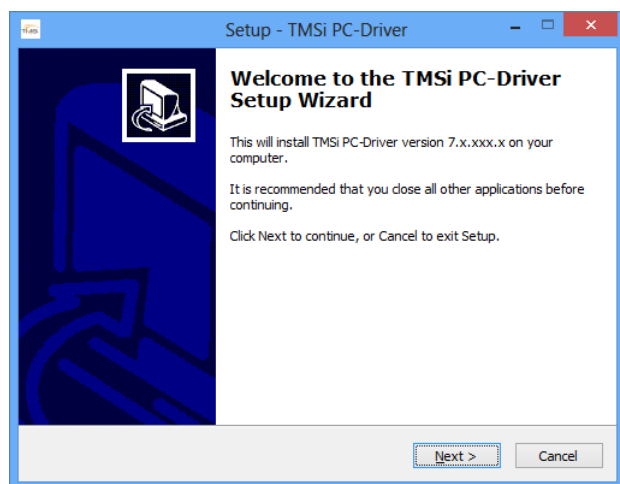


Step 4: Accept the License Agreement and click ‘Next’.

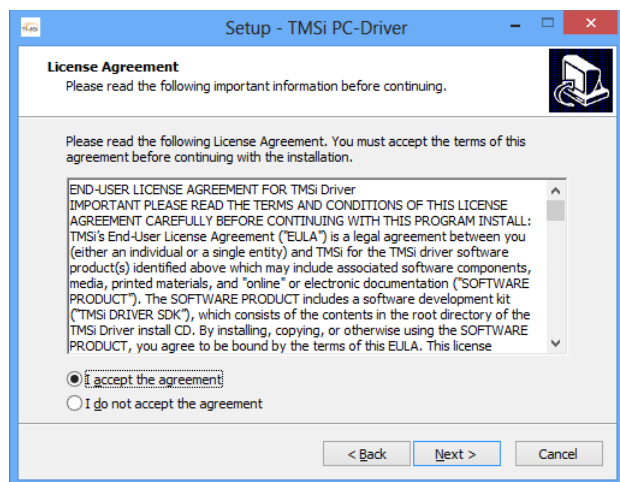

Step 5: Select Installation Folder and click 'Next'.

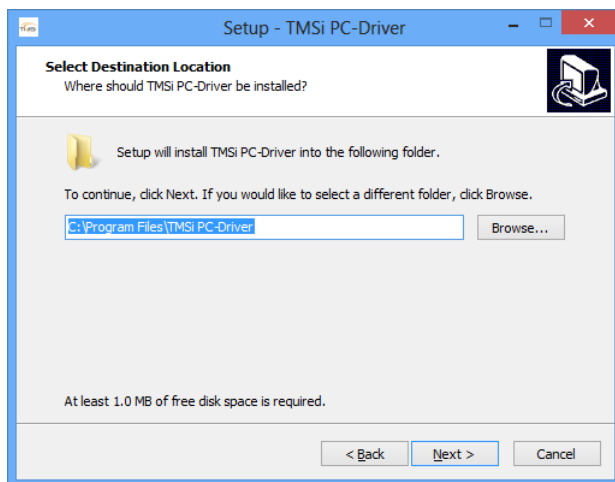


Step 6: 'Ready to install', click 'Install'.

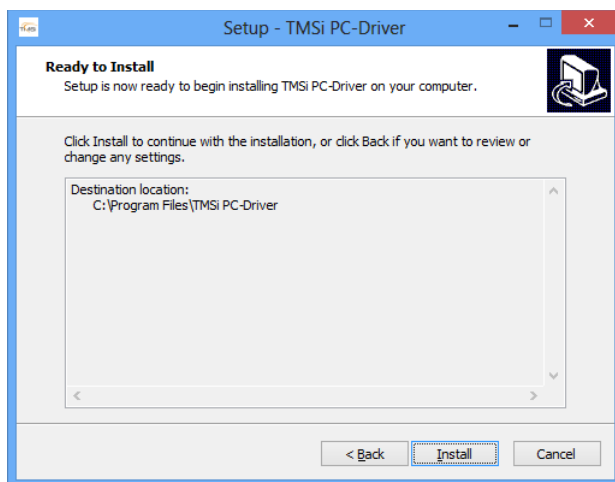

Step 7: 'Device Driver Installation Wizard' pops up, click 'Next'.



Step 8: If 'Windows Security' pops up, then tick the box 'Always trust software from "TMS INTERNATIONAL B.V."' and click 'Install'.

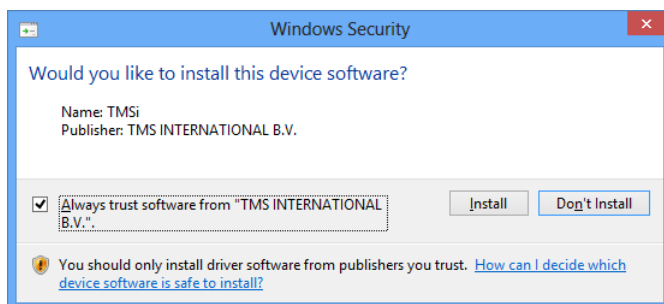


Step 9: 'Completing the Device Driver Installation Wizard', click 'Finish'.



Now the Microsoft Visual C++ Redistributable will be installed, this requires no user action.

Step 10: Setup is complete, click 'Finish'.

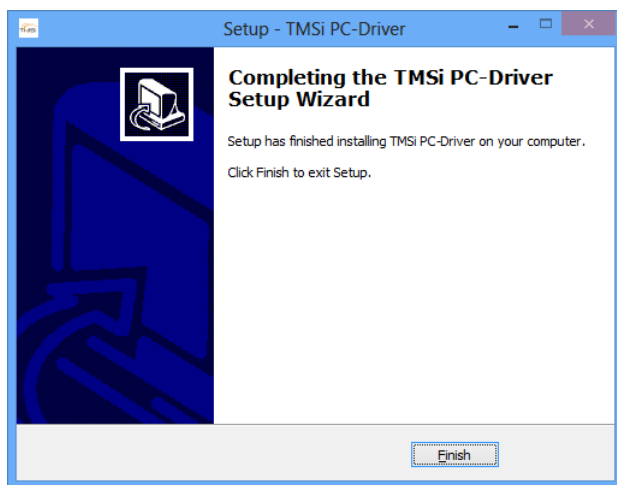

Thank you for installing the TMSi PC-driver! If you have any questions do not hesitate to contact us at support@tmsi.com.

Copyright © 2014 TMSi. All rights reserved.